







- Simple preparation; Dilute it 10 times
 with distilled water
- Methanol-Free
 - → Easy waste liquid treatment
- Fast electrophoresis transfer
- Available for Wide range of molecular weight

Fast electrophoresis transfer

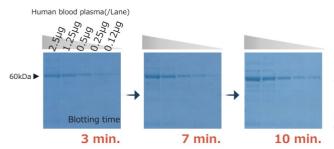

- ightarrow Fast blotting within 10 min is possible.
- - membrane: P + memebranes (Pore size:0.2µm)

XX Available for wide range of molecular weight (~200kDa)

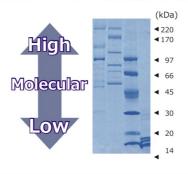

- $\rightarrow \textit{EzFastBlot}$ is available for wide range of molecular weight, from low to high.
 - $\ensuremath{\ensuremath{\%}}$ Each data is membrane after blotting and dyeing with CBB stain.
 - membrane: P + memebranes (Pore size:0.2µm)

Product	AE-1465 EzFastBlot
Code	2332590
Package	1 Bottle
Volume	500 mL (10x concentrated solution)
Use	For blotting x 50 times
Main component	Tris, Amino acids
Storage	at room temperature for 1 year

